

Alpha-Theta Synchrony Neurofeedback Module for Nexus

2 channels

Sensor "DC Offset"

Asymmetry Monitors

Raw Electrocardiogram

Heart Rate Variability
Power Spectrum
(Autonomic balance)


Pulse

Heart Rate Variability

HRV - Low Freq. %

Temp

GSR

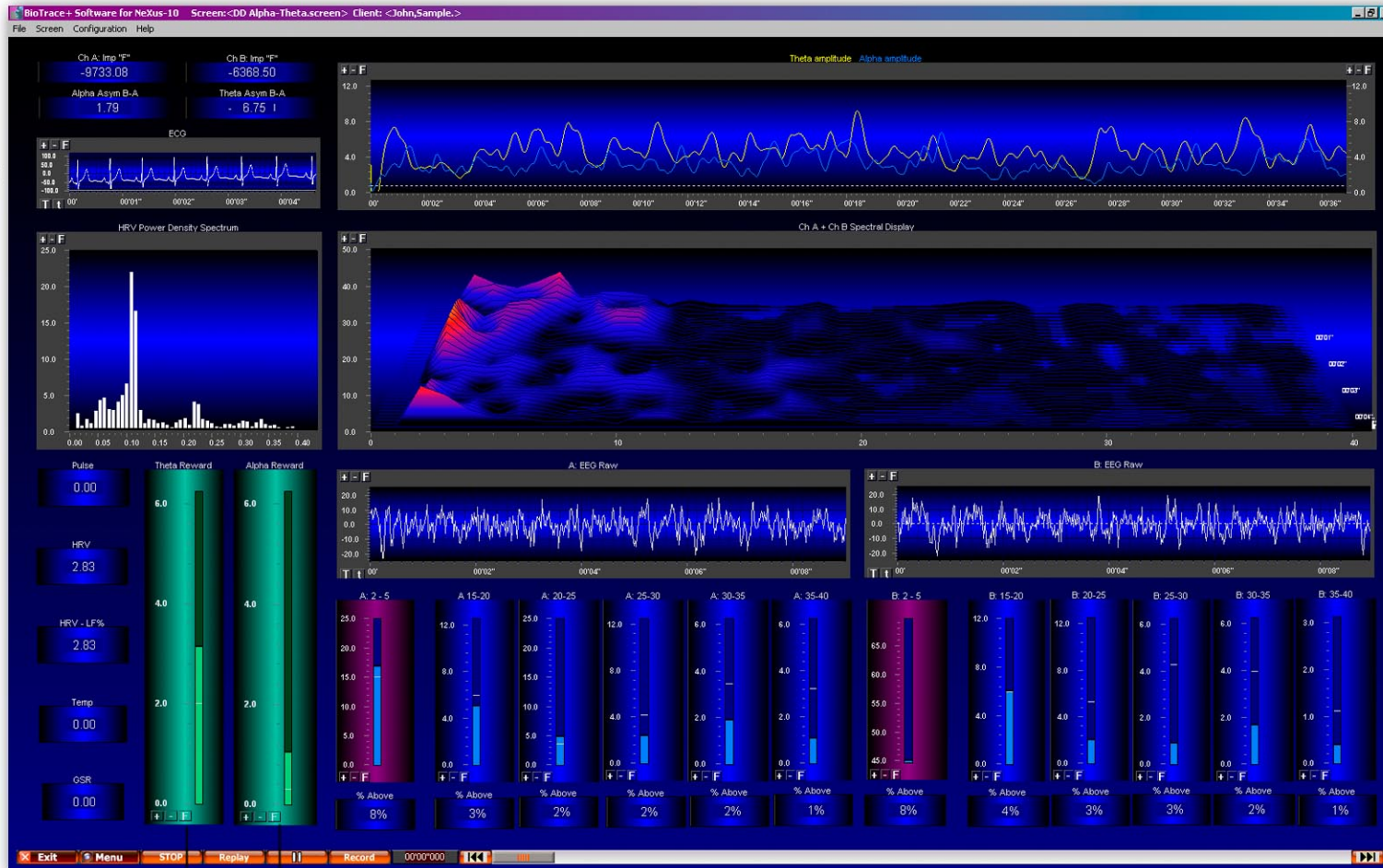

Alpha-Theta
Crossover

Spectral Display
A + B

Raw EEG

6 Programmable
Inhibit bins per
Channel

% Time Above
Threshold

Completely programmable alpha reward. Comes preset to reward a 3 Hz bin centered on 10 Hz.

Completely programmable theta reward. Comes preset to reward a 3 Hz bin from 5 Hz to 8 Hz.

Comes with : Nine theta feedback sounds, namely, 5, 6 & 7 Hz binaural beats centered on 40, 80 & 120 Hz. Comes with several choices for alpha feedback, namely, plain "pink" noise, pink noise embedded with alpha, etc. You can also use your own mp3 or wave files. Low frequency inhibit uses Nexus "Birds" sound. High frequency inhibits are set to a "blip" sound to help control rogue excursions. Registered users are entitled to free technical support. Upgrades are always free.

Custom Nexus applications can be handcrafted to your specifications. For more information call "A Matter of Mind ®" at (408) 984-3333 or email: mind@growing.com